




DICOM Conversion and Integration on touchscreen for medical use


Digital recording system integrating images and motion video of medical records on touchscreen




Dermatology

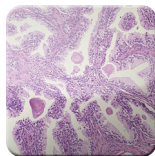


Pathology

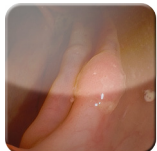

Endoscopy

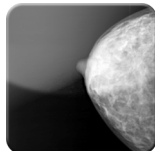

Mammography

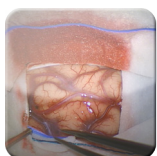


Neurosurgery



Ophthalmology



Ultrasounds

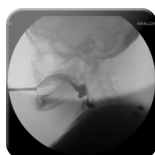

Fluoroscopy

Easy to use and intuitive

ETIAM DICOM Izer is the most ergonomic and complete solution for capturing, edit images and video sequences in your specialty, and to distribute them in DICOM format.

Extended customizability

You can easily fit ETIAM DICOM Izer to your own specialty workflow and thus automate its use and finely integrate it into your imaging network.

Compliance with standards


In compliance with the DICOM standard, IHE recommendations and the HIPAA regulation, ETIAM DICOM Izer enables to secure patient data and enhance traceability.

Increased productivity at a lower cost

ETIAM DICOM Izer reduces costs since it integrates directly into the existing infrastructure (PACS and patient record). Since you access images and video directly through the patient record, you save time. Also, with ETIAM DICOM Izer you do not have to change your work habits. Then, you have more time to focus on your patients.

Recommended Technical Specifications


- Touch screen: 15", 1024 x 768, 400 cd/square cm
- Video inputs: Full resolution PAL, NTSC (S-Video, Composite)
- Audio: Jack 3.5 mm
- Processor: Pentium M 1.6 GHz; RAM: 1 GB
- Storage: CD/DVD, 80GB internal disk
- Dimensions: W: 398 x H: 330 x D: 93 mm
- Weight: 8 Kg
- Mounting: VESA 75
- Peripherals: Remote control (foot pedal, VTR)

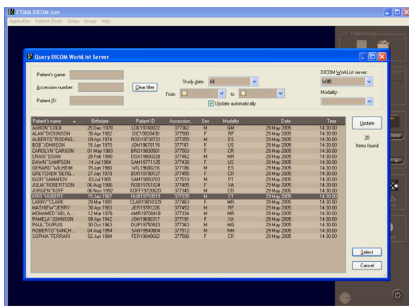

Select a patient



Capture from the endoscope

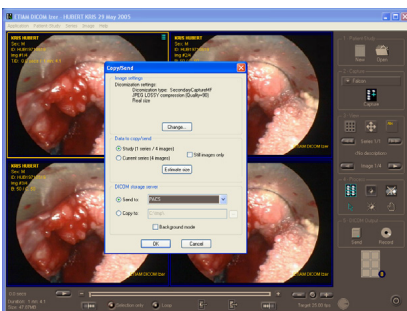

Distribute

Features

- Touch screen interface
- Multiple switchable video inputs
- Remote control (endoscope triggers, foot pedal, etc.)
- Standalone (Print, CD/DVD, USB key) and/or connected
- (DICOM Worklist, Store to PACS, etc.) solution
- Record video, audio and still images simultaneously

Video edition (video trimming)

Multiple output formats provided (from AVI, MPEG to

DICOM formats, including Multiframe and MPEG)

Certification/markings

CE/FCC Class B, UL 60601-1, EN 60601-1

Supported Endoscopes (partial)

- FUJINON
- OLYMPUS
- PENTAX
- SOPRO-COMEG
- STORZ
- Stryker
- WOLF
- etc.

Benefits

Security and quality of care: Create a link between moving images and all other patient clinical information

Extend your medical report content from still images only to video (for some studies digital video recording is a demand)

Distribute moving images to referring doctors on a CD or DVD with a powerful and intuitive DICOM viewer

ETIAM Corp.

100 Main Street
Suite 10
Maynard, MA 01754
USA

Toll Free: (877) 384-2662
Tel./Fax: (617) 395-5809
E-mail: sales-us@etiam.com

www.etiam.com

ETIAM Group

La Palmeraie
ZA de la Hallerais
11 rue du Bois de Soeuvers
35770 Vern-sur-Seiche
France

Tel.: +33 (0)2 99 14 33 88
Fax: +33 (0)2 99 14 33 80
E-mail: sales-fr@etiam.com

ETIAM Asia

Bld 6, No. 1279, Zhongshan
West Road
200051 Shanghai
Office 251
China

Tel.: +86 0215-116-8268
Fax: +86 1561-830-8300
E-mail: sales-cn@etiam.com

More information on

www.etiam.com

This document is not contractual and this offer is subject to changes without prior notice.
All trademarks and registered trademarks herein are the property of their respective owners.